


Type	CPAP	Auto-CPAP	Auto Bi-level
Model	C-20C	C-20A	B-25A
Respiratory mode	CPAP	CPAP, APAP	CPAP, S, AutoB
<b>Parameters</b>			
Pressure	4-20cmH <sub>2</sub> O	4-20cmH <sub>2</sub> O	4-20cmH <sub>2</sub> O
Islop	-	-	1-6 Level
95th Pressure titration	-	Available	Available
Max PS	-	-	3-8cmH <sub>2</sub> O
Max Press	-	4-20cmH <sub>2</sub> O	-
Min Press	-	4-20cmH <sub>2</sub> O	-
IPAP	-	-	4-25cmH <sub>2</sub> O
EPAP	-	-	4-23cmH <sub>2</sub> O
Max IPAP	-	-	4-25cmH <sub>2</sub> O
Min EPAP	-	-	4-23cmH <sub>2</sub> O
Auto-On	On/Off	On/Off	On/Off
Auto-Off	On/Off	On/Off	On/Off
Humidifier	0-5levs	0-5levs	0-5levs
Auto-air compensation	Available	Available	Available
Auto-adjustment altitude	Available	Available	Available
Noise Level	Less than 30dBA	Less than 30dBA	Less than 30dBA
<b>Configuration</b>			
Screen	3.5"Colorful LCD	3.5"Colorful LCD	3.5"Colorful LCD
SpO2 Kit	Optional	Optional	Optional
Data Storage	TF Card	TF Card	TF Card
Dimension	280*200*112mm	280*200*112mm	280*200*112mm
Weight	2.2Kg	2.2Kg	2.2Kg




**ResPlus**  
CPAP/APAP series  
For Sleep Apnea Therapy

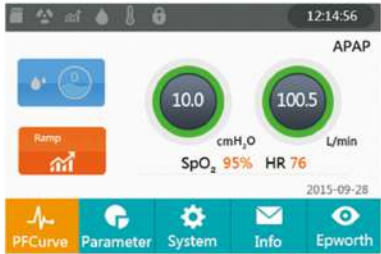

### 3.5" Color LCD Screen

- All essential information is easy to read
- User-friendly interface enables fast setup & easy control

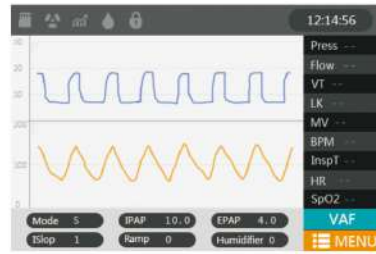

### Real-time Data Monitoring

- Real-time key data monitoring ensures an effective therapy

### SpO2 Monitoring

- SpO2 Kits(optional) enable real-time SpO2 monitoring


### Ultra-Quiet Design

- High-end Swiss blower guarantees a low dBA (< 30)



### Safety Design with Heating Features

- 5 temperature levels available for humidifier
- Back flow prevention design

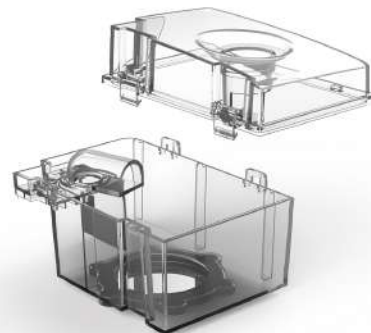

### Detachable Design

- Detachable humidifier is easy to carry, clean and fit into different settings



## ResPlus C-20C MODE: CPAP

### Pressure Relief Technology

Reduce the exhalation pressure to maximize therapy comfort

### Air Leak Compensation

Automatic leakage detection and air compensation

### Pressure Compensation

Automatic compensation for different pressure scenarios

## ResPlus C-20A MODE: APAP, CPAP

### Accurate Detection

High-end sensors with innovative algorithms to identify various respiratory scenarios precisely and respond quickly

### Auto-adjusting Pressure

Automatically adjust the pressure to optimal level in accordance with different respiratory events

### Pressure Titration

Recommend the best therapy pressure setting based on patient's actual data

### Real Time Alerts

Real time audible alerts for power failure, air leak, humifailure, etc

## ResPlus B-25A MODE: CPAP, S, AutoB (Customized & Intelligent)

### Auto Bilevel

Auto-adjustment of IPAP and EPAP in accordance with patient's respiratory events

### Auto-sensibility

Isens and Esens are automatically adopted to individual's treatment

